

HoverMatt® Air Transfer System

From HoverTech International

User Manual



T: 1300 106 003

F: 07 3036 6625

E: info@hmsmedical.com.au

www.hmsmedical.com.au


HOVERMATT®
www.HoverMatt.com

HOVERMATT® USER MANUAL

Table of Contents

Symbol References	02
Intended Use and Precautions	03
HOVERMATT® Introduction	04
Part Identification - Mattress	05
Part Identification - Air Supply	06
Air Supply Keypad Functions	07
Instructions for Use	08
Product Specifications/Required Accessories	09, 10
Electromagnetic Compatibility Chart	11 - 14
Cleaning	15
Preventative Maintenance/Infection Control	16
Frequently Asked Questions	17
AIR SUPPLY	18
Part Identification	19
Power Cord / Clamp Replacement	20
Handle Replacement	21
Feet or Bumper Replacement	22
Hose Removal	23
Air Filter and Air Filter Cover Replacement	24
Dust Cover/Hose Attachment Snap Replacement	25
Metal Cover Replacement	26
Cord Strap Replacement	27
Troubleshooting	28
General System Information	
Components Parts List	29
Warranty Statement	30, 31
Returns and Repairs	32

HOVERMATT® USER MANUAL

Symbol Reference




Attention! Please read accompanying documents.



This End Up



Type BF Applied Part

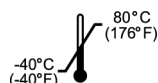

Temperature



Declaration of Conformity to Medical Device Directive



Humidity



Functional Earth (Ground)



Date of Manufacture



Alternating Current



Keep Dry



Underwriters Laboratory Agency Approval

120 V~:

Medical Equipment with respect to electrical shock, fire and mechanical hazards only in accordance with UL 60601-1, IEC/EN 60601-1, CAN/CSA C22.2 No. 601.1

230 V~:

Medical Equipment with respect to electrical shock, fire and mechanical hazards only in accordance with UL 60601-1, IEC 60601-1-2 CAN/CSA C22.2 No. 601.1

HOVERMATT® USER MANUAL

Intended Use and Precautions

- **Indications:**
Patients unable to assist in their own lateral transfer. Patients whose weight or girth poses a potential health risk for the caregivers responsible for repositioning or laterally transferring said patients.
- **Contraindications:**
Patients who are experiencing thoracic, cervical or lumbar fractures that are deemed unstable, unless using in conjunction with a spinal board on top of the HoverMatt mattress.
- **Intended Care Settings:**
Hospitals, long term or extended care facilities
- **Precautions:**
Caregivers must verify that all caster brakes have been engaged prior to transfer.
Additional caregivers are recommended when moving a patient over 750 lbs / 340kg.
Never attempt to move a patient on an uninflated HoverMatt mattress.
Route the power cord in a manner to ensure freedom from hazard.
Avoid blocking the air intakes of the Air Supply.
Never leave patient unattended on an inflated device.

Use this product only for its intended purpose as described in this manual. Only use attachments and / or accessories that are authorized by HoverTech International.

When transferring to a low air loss bed, set the bed mattress air flow to the highest level for a firm transfer surface.

When using the HoverMatt in the MRI environment, a 25 ft. specialty MRI hose is required (available for purchase).

WARNING: For safety, always use two people during patient transfer.

CAUTION: Avoid electric shock. Do not open Air Supply.

WARNING: Reference product specific user manuals for additional operating instructions.

HoverJack® Air Patient Lift system is not UL classified

Introduction

Using the HoverMatt® Air Transfer System

The HoverMatt® air transfer system is used to assist with lateral transfers and repositioning anywhere in the hospital. It is radiolucent and artifact free, so patients may remain on the HoverMatt® transfer mattress for all ancillary procedures, including Diagnostic Radiology, Nuclear Medicine, MRI, or Radiation Therapy. The patient weight limit for the HoverMatt® air transfer system is 1,200lbs / 544kg. It is available in four different widths to accommodate the body mass of the patient.

The Principle of the HoverMatt® Technology


After the patient is placed on the HoverMatt® transfer mattress, low pressure air from the small Air Supply will inflate the mattress. At the same time the air is supporting the patient, the air is escaping from the perforations in the underside of the HoverMatt® transfer mattress. The escaping air acts as a lubricant to reduce friction, which facilitates effortless transfers. With less force needed to transfer a patient, there is less physical effort and strain expended by caregivers, which results in a reduction of workers' compensation injuries. The HoverMatt® air patient transfer system requires the caregiver to exert a force of approximately 15% of the patient's body weight for the transfer.

The Purpose of HoverMatt® Technology

Consistent utilization of the HoverMatt® air transfer system dramatically reduces back injuries to staff that are caused by lateral transfers and repositioning. In addition, fewer staff members are required to perform these tasks and a more comfortable transfer is provided for the patient.


HOVERMATT® USER MANUAL

Part Identification - Mattress

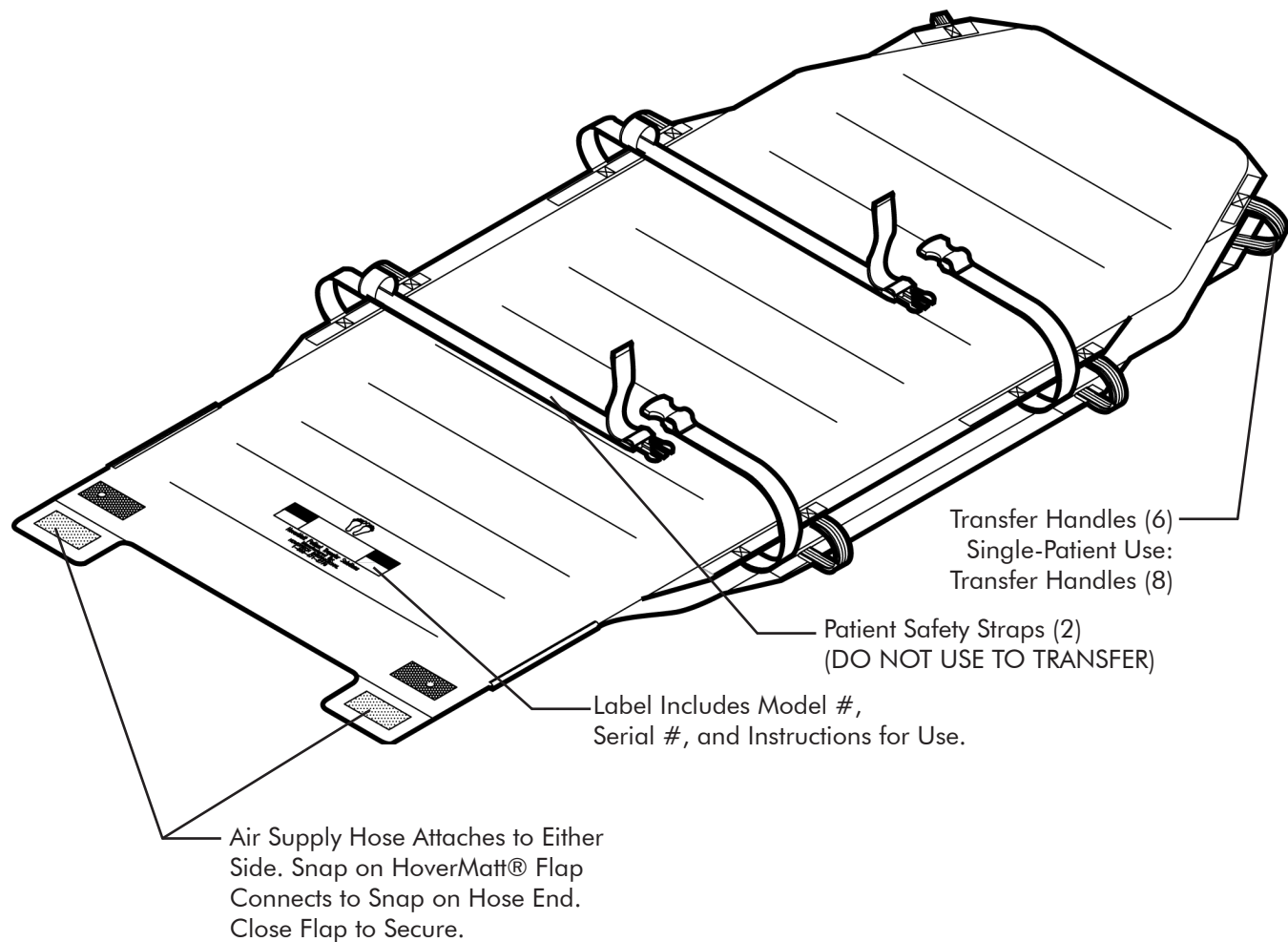

HOVERMATT® USER MANUAL



Part Identification - Air Supply

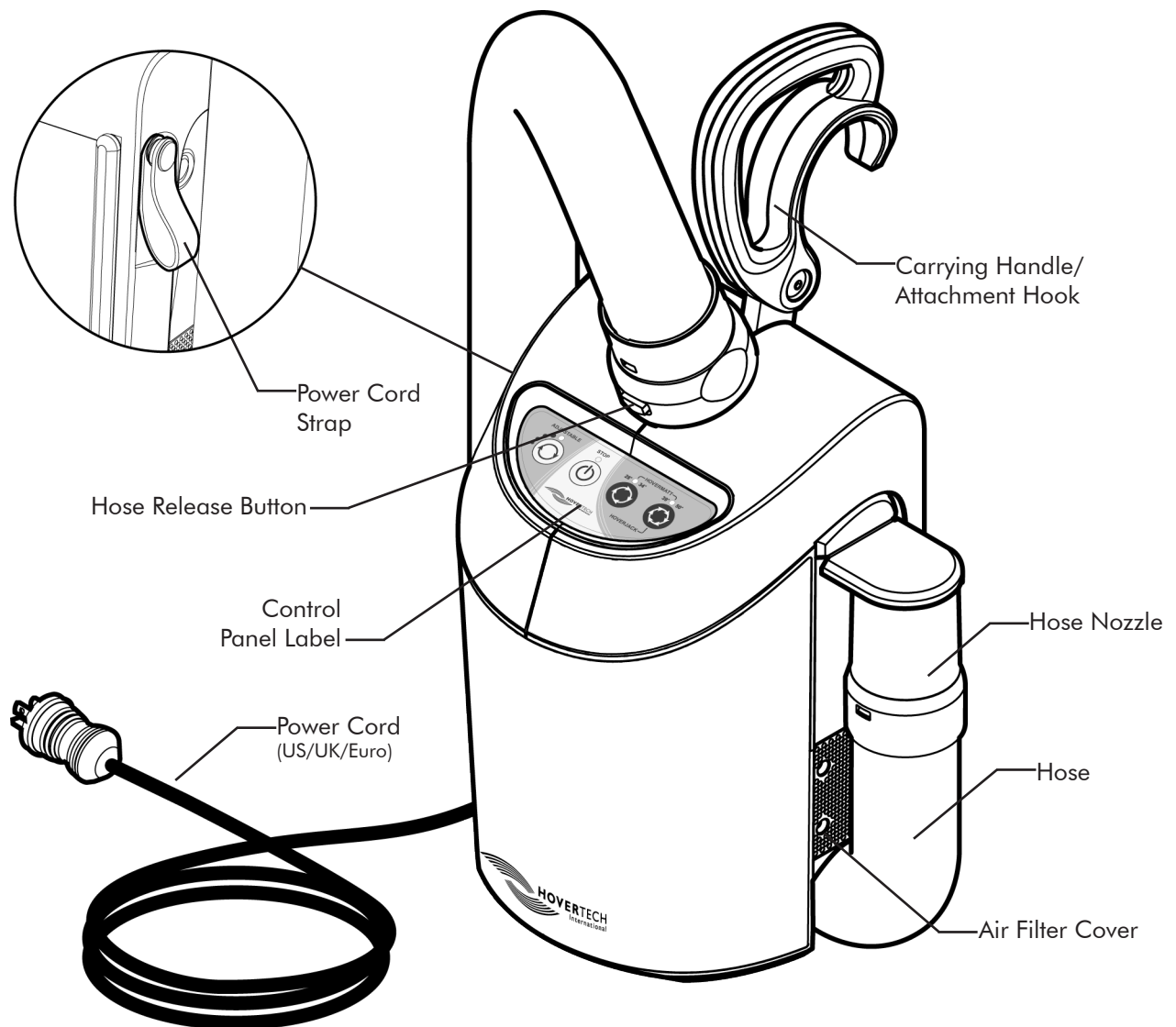

WARNING: The HTAIR is not compatible with DC power supplies. The HTAIR is not for use with the HoverJack Battery Cart.

HOVERMATT® USER MANUAL

Air Supply Keypad Functions

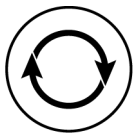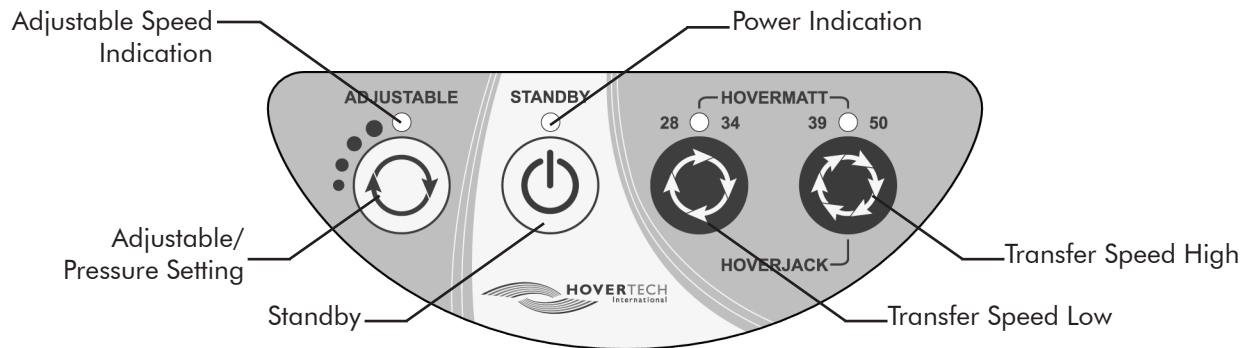


The Adjustable Keypad function has four different settings. Pressing the button once will result in the lowest inflation setting available. A second press of the button increases the air pressure and rate of inflation. Pressing the button a third time will again increase the rate of inflation. A fourth press of the button results in the highest inflation rate and air pressure available for HoverTech Accessories. The STANDBY button may be pressed at any time to cease all air flow.

NOTE: The LED will indicate the inflation speed by the number of flashes (i.e. two flashes equals the second inflation speed).

All of the settings in the Adjustable range are substantially lower than the HoverMatt® and HoverJack® settings. The Adjustable function is not to be used for transferring; it is only for use with HoverTech Accessories, which require a lower pressure for slow inflation.



Standby: Used to stop inflation/air flow.

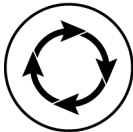

HoverMatt® 28 / 34: For use with 28" & 34" HoverMatt® Air Transfer Mattresses.



HoverMatt® 39 / 50 & HoverJack®: For use with 39" & 50" HoverMatt® Air Transfer Mattresses and 32" & 39" HoverJack® Air Patient Lifts.


HOVERMATT® USER MANUAL

Instructions for Use

1. Patient should be in a horizontal position for transfer/repositioning on the appropriate width mattress. [28"w (71 cm) – 34"(86 cm) – 39"w (99 cm) – 50"w (127 cm)]

Regardless of ease of patient movement, for safety, always use a minimum of two caregivers for the transfer.

2. Place the HoverMatt® transfer mattress underneath patient using log-rolling technique and attach patient safety straps. Whatever the patient is lying on to keep the bed mattress clean can be placed on top of the HoverMatt® transfer mattress to help keep it clean.
3. Plug the Air Supply power cord into an electrical outlet.
4. Insert the Air Supply hose nozzle into mattress in one of the two entries located parallel to the foot end, and snap in place.
5. Be sure transfer surfaces are as close as possible and brake the wheels.
6. If possible, transfer from a higher surface to a lower surface.
7. Turn on Air Supply by choosing the appropriate speed based on the mattress width.

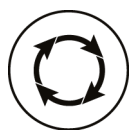

HOVERMATT®

Use for sizes: 28 & 34



HOVERMATT®

Use for sizes: 39 & 50

8. Grasp transfer handles and pull patient on an angle, either head first or feet first, until patient is in desired position.
9. Ensure that the patient is centered on the receiving equipment prior to deflation, especially if the width of the equipment receiving the patient is less than the width of the transfer mattress.
10. Press Standby to deflate mattress and employ the bed/stretchers rails.

NEVER LEAVE PATIENT UNATTENDED ON AN INFLATED HOVERMATT® AIR TRANSFER SYSTEM.

Product Specifications/Required Accessories

- Classification:
Not for use in the presence of flammable anesthetics or in a hyperbaric chamber or oxygen tent.
Type of Protection against electric shock: CLASS I EQUIPMENT
Degree of protection against electric shock: TYPE BF APPLIED PART
Protection against ingress of water: Ordinary (not protected).
Mode of operation: CONTINUOUS OPERATION
To remove supply mains, unplug equipment from wall.
- Patient Weight Limit: 1,200lbs / 544kg
- Use Temperature: 50° to 104° F (10° to 40° C)
- Use Humidity: 10% to 70% Non-Condensing
- Storage/Shipping Temperature: -40° to 176° F (-40° to 80° C)
- Storage/Shipping Humidity: 10% to 70% Non-Condensing
- Power Input: 120 V~, 60 Hz, 10 A (North American version)
230 V~, 50 Hz, 6 A (European Version)
- Air Supply Dimensions: 12.5 x 7 x 7 inches (31.75 x 17.8 x 17.8 cm)
- Air Supply Weight: 11 lbs. (5 kg)
- Air Supply Material: Fire Retardant ABS/Stainless Steel
- Power Cord Length: 15 feet (457 cm)

HOVERMATT® USER MANUAL

Product Specifications/Required Accessories

	RE-USABLE HOVERMATT® AIR TRANSFER SYSTEM	SINGLE-PATIENT USE HOVERMATT® AIR TRANSFER SYSTEM
Material:	Nylon Twill	Nylon Twill & Non-Woven
Construction:	RF-Welded	Sewn
Width:	28" (71 cm), 34" (86 cm), 39" (99 cm) or 50" (127 cm)	34" (86 cm), 39" (99 cm) or 50" (127 cm)
Length:	78" (198 cm)	78" (198 cm)

Required Accessory:

HoverTech International Air Supply

Part# HTAIR1200 (North American Version)

Part# HTAIR2300 (European Version)

All HoverTech International Products are Latex-Free.

For a full product listing go to www.HoverMatt.com.

Electromagnetic Compatibility Chart

For HTAIR-2300 ONLY

Guidance and Manufacturer's Declaration – Electromagnetic Emissions		
<p>The HoverTech International Air Supply is intended for use in the electromagnetic environment specified below. The customer or the user of the HoverTech International Air Supply should ensure that it is used in such an environment.</p>		
Emissions Test	Compliance	Electromagnetic Environment Guidance
RF emissions CISPR11	Group 1	The HoverTech International Air Supply uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR11	Class A	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC-61000-3-3	Complies	The HoverTech International Air Supply is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.


Electromagnetic Compatibility Chart

For HTAIR-2300 ONLY

Guidance and Manufacturer's Declaration – Electromagnetic Immunity			
The HoverTech International Air Supply is intended for use in the electromagnetic environment specified below. The customer or the user of the HoverTech International Air Supply should ensure that it is used in such an environment.			
Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment-Guidance
Electrostatic Discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast Transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input-output lines	± 2 kV for supply mains ± 1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	< 5% U_T (> 95% dip in U_T) For 0,5 cycle 40% U_T (60% dip in U_T) For 5 cycles 70% U_T (30% dip in U_T) For 25 cycles < 5% U_T (> 95% dip in U_T) for 5 seconds	< 5% U_T (> 95% dip in U_T) For 0,5 cycle 40% U_T (60% dip in U_T) For 5 cycles 70% U_T (30% dip in U_T) For 25 cycles < 5% U_T (> 95% dip in U_T) for 5 seconds	Mains power quality should be that of a typical commercial or hospital environment. If the user of the HoverTech International Air Supply requires continued operation during mains interruptions, it is recommended that the HoverTech International Air Supply be powered from an uninterruptible power supply or a battery.
Power Frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical commercial or hospital environment.
NOTE: U_T is the AC mains voltage prior to application of the test level			

Electromagnetic Compatibility Chart

For HTAIR-2300 ONLY

Guidance and Manufacturer's Declaration –Electromagnetic Immunity			
The HoverTech International Air Supply is intended for use in the electromagnetic environment specified below. The customer or the user of the HoverTech International Air Supply should ensure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 V	Portable and mobile RF communications equipment should be used no closer to any part of the HoverTech International Air Supply, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	<p>Recommended separation distance $d = 1.2 \sqrt{P}$</p> <p>$d = 1.2 \sqrt{P}$ 80 to 800 MHz</p> <p>$d = 2.3 \sqrt{P}$ 800 MHz to 2.5 GHz</p> <p>Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b.</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
<p>NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.</p> <p>NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.</p> <p>a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, and electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the device.</p> <p>b Over the frequency range 150 kHz to 80 MHz, the field strengths should be less than 3 V/m.</p>			

Electromagnetic Compatibility Chart

For HTAIR-2300 ONLY

Recommended separation distances between portable and mobile RF communications equipment and the HoverTech International Air Supply			
<p>The HoverTech International Air Supply is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the HoverTech International Air Supply can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the HoverTech International Air Supply as recommended below, according to the maximum output power of the communications equipment.</p>			
Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	15 kHz to 80 MHz $d=[3.5/V_1]\sqrt{P}$	80 MHz to 800 MHz $d=[3.5/V_1]\sqrt{P}$	800 MHz to 2.5 GHz $d=[7/E_1]\sqrt{P}$
0.01	0.12	0.12	0.23
0.10	.38	.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23
<p>For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.</p> <p>NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.</p> <p>NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.</p>			

Cleaning

Reusable HoverMatt® Transfer Mattress

The reusable HoverMatt® transfer mattresses are constructed of nylon twill.

In between patient use, the HoverMatt should be wiped down with a 100:1 bleach solution (100 parts water: one part bleach) or cleaning solutions your hospital uses for medical equipment disinfections.

If the HoverMatt becomes badly soiled, it can be laundered in a washing machine with a 160° degrees Fahrenheit (65° C) maximum water temperature. A 100:1 bleach solution may be used (100 parts water: one part bleach) during the wash cycle.

The mattress should be air dried if possible. Air drying can be expedited by using the Air Supply unit to circulate air through the inside of the mattress. If using a dryer, the temperature setting should be set on the coolest setting. Drying temperature should never exceed 115° degrees Fahrenheit (46° C). The backing of the nylon is polyurethane, and will begin to deteriorate after repeated high temperature drying. The Double-Coated HoverMatt® transfer mattress should not be put in the dryer.

To help keep the HoverMatt® transfer mattress clean, HoverTech International recommends the use of their disposable or reusable protector sheets. A bed sheet may also be used.

Single-Patient Use HoverMatt® Transfer Mattress

The Single-Patient Use HoverMatt® transfer mattress is not intended to be reprocessed. If the Single-Patient Use HoverMatt® transfer mattress is laundered, the perimeter seam thread will dissolve.

Air Supply Cleaning and Maintenance

In between patient uses, the Air Supply can be cleaned by wiping down using a damp cloth with soap and water or mild neutral detergent. Dry using a clean, dry cloth or disposable paper towel.

*Do not spray cleaners or liquids directly on the Air Supply.

NOTE: CHECK YOUR LOCAL/STATE/FEDERAL/INTERNATIONAL GUIDELINES BEFORE DISPOSAL.

Preventive Maintenance

Prior to use, a visual inspection should be performed on the HoverMatt® air transfer system to insure the Air Supply power cord is not frayed or nicked, and that there is no visual damage that would render the Air Supply unusable. The HoverMatt® transfer mattress should have all of its safety straps and handles (reference Page 4 of the manual for all appropriate parts). There should be no tears or holes that would prevent the HoverMatt® transfer mattress from inflating. If any damage is found that would cause the system not to function as intended, the HoverMatt® air transfer system should be removed from use and returned to HoverTech International for repair (see Page 32 Repairs & Returns).

The Air Supply has air filters on either side of the motor. These filters can be accessed by removing the small screws holding the filter cover in place. Filters should be cleaned by holding under warm running water. Allow to air dry. As preventive maintenance, filter cleaning should be performed monthly.

Infection Control

HoverTech International offers superior infection control with our heat-sealed reusable HoverMatt® transfer mattress. This unique construction eliminates the needle holes of a sewn mattress which can be potential bacterial entry ways. Additionally, the heat-sealed, double-coated HoverMatt® transfer mattress offers a stain and fluid proof surface for easy cleaning. A Single-Patient Use HoverMatt® transfer mattress is also available.

Whatever the patient is lying on to keep the hospital bed clean may be placed on top of the HoverMatt® transfer mattress to help keep it clean. If desired, the protector sheet or disposable sheet may be used to cover the transfer mattress (available for separate purchase). This sheet may also be placed under the transfer mattress, when log rolling the patient, to prevent the bottom of the transfer mattress from coming in contact with the hospital bed sheets.

If the HoverMatt® transfer mattress is used on an isolation patient, the hospital should employ the same protocols/procedures it utilizes for the bed mattress and/or for the linen in that patient room.

Frequently Asked Questions

1. *What is the weight limit of the HoverMatt® air transfer system?*

The patient weight limit for the HoverMatt® air transfer system is 1,200lbs / 544kg. It is available in four different widths to accommodate the body mass of the patient.


2. *Is the HoverMatt® air transfer system latex free?*

Yes.

3. *Is the HoverMatt® air transfer system FDA Approved?*

HoverTech International is registered with the FDA. The HoverMatt® air transfer system is listed as a Class II device and is exempt from pre-market notification and approval.

AIR SUPPLY

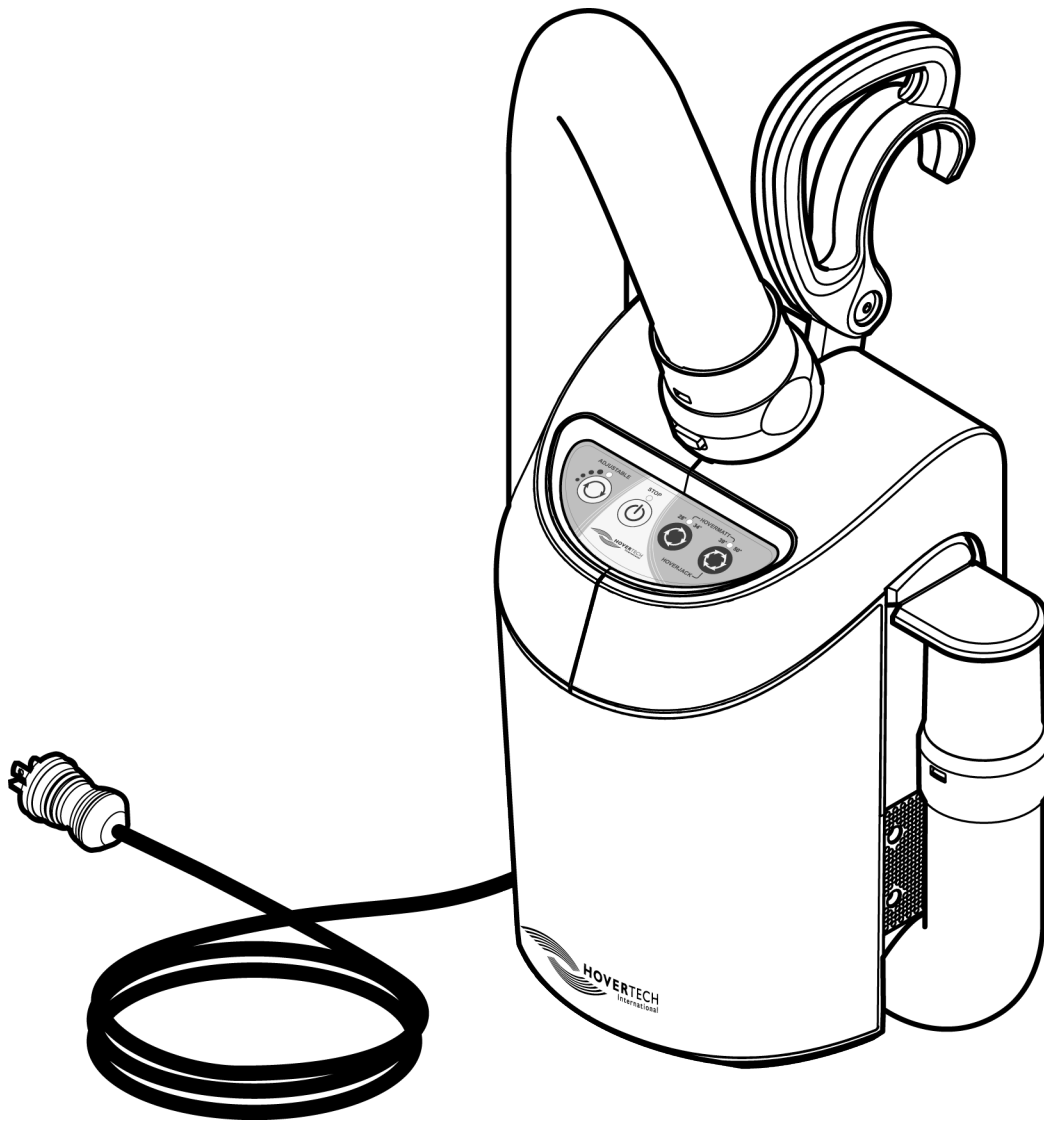


NO USER SERVICEABLE PARTS.

Only qualified service personnel shall perform repairs on the HoverTech International Air Supply.


HOVERMATT® USER MANUAL

Part Identification

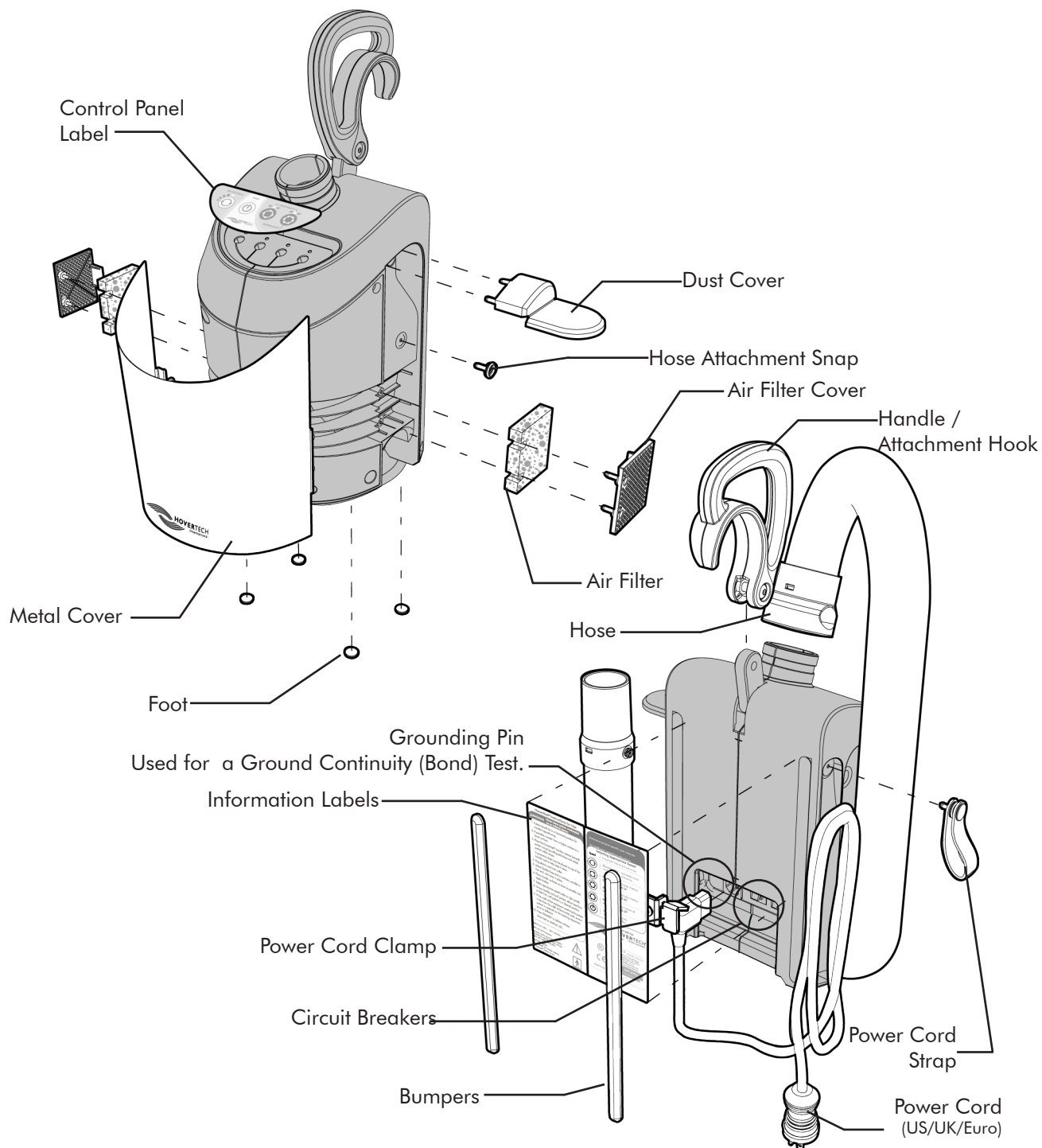


HOVERMATT® USER MANUAL

Power Cord / Clamp Replacement

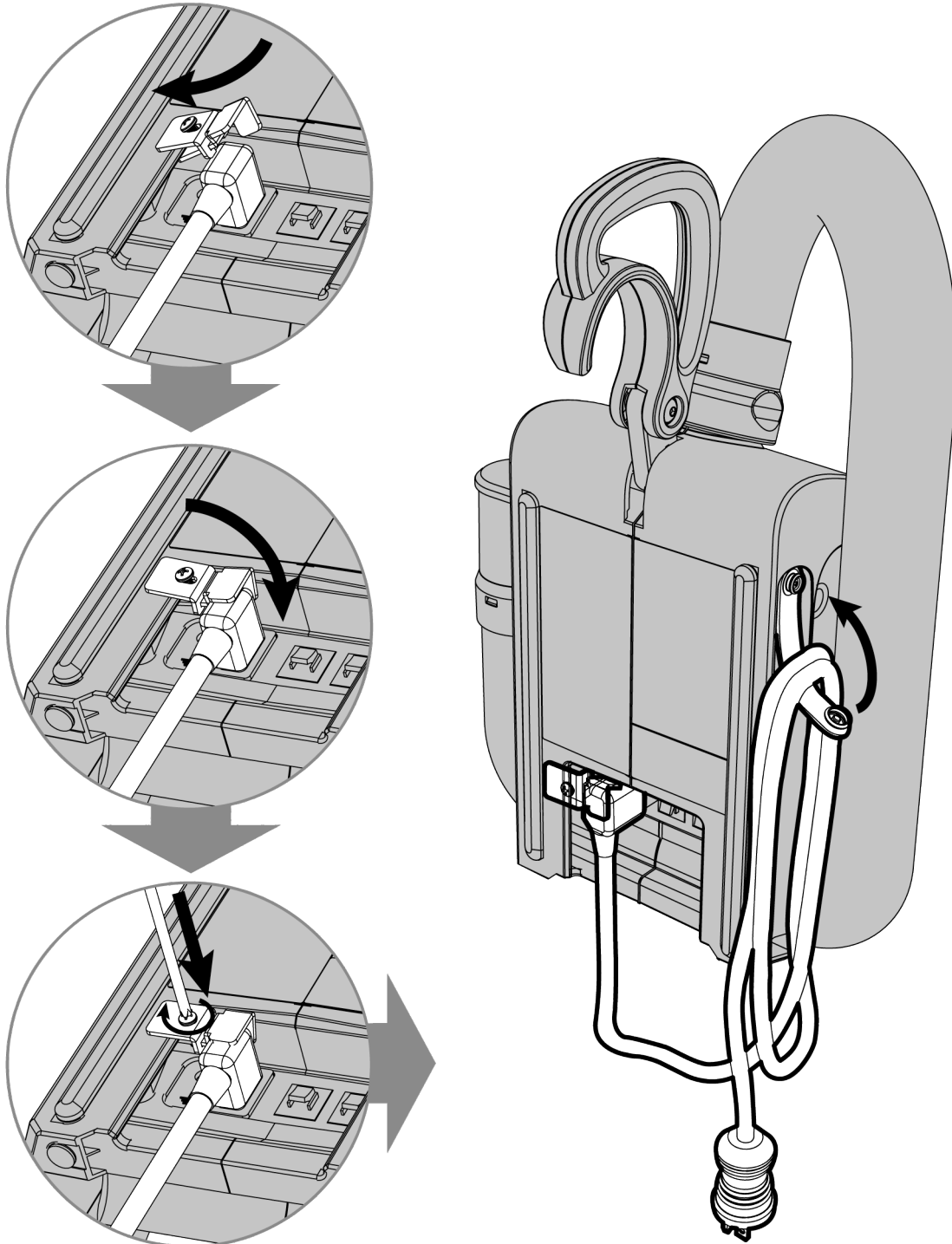


Handle Replacement

1. Remove the damaged handle by unscrewing the socket head screw from the barrel nut using two 5/32" allen wrenches as shown.
2. Attach the new handle by reversing the process. When tightening the screw be sure that the handle can rotate easily. The screw is treated with thread lock to secure it in place.

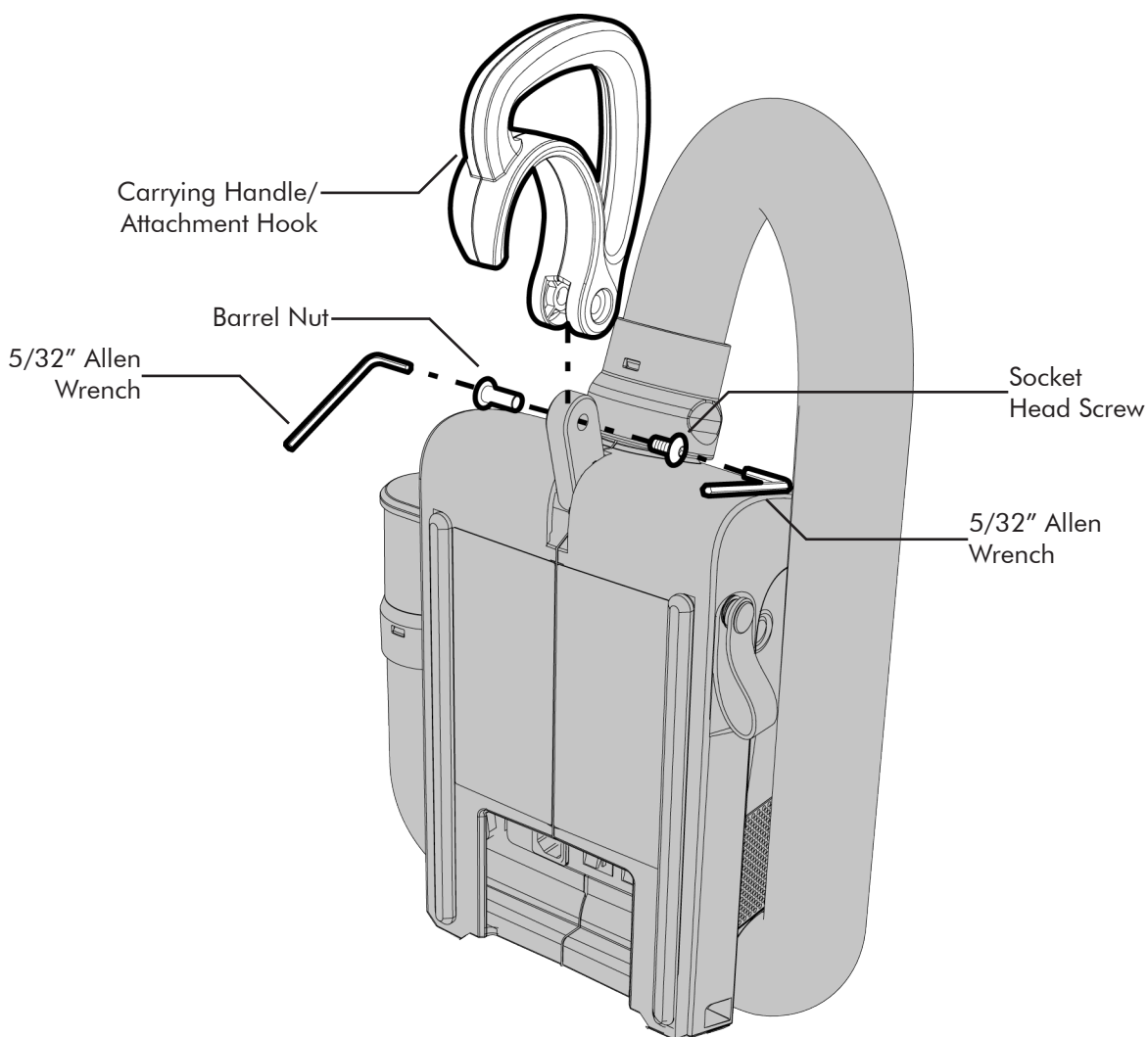


Feet or Bumper Replacement

1. The feet and bumpers are held in place by a self-adhesive backing. Use a small, flat bladed screwdriver to pry up an edge and gently remove the foot or bumper.
2. Clean surface thoroughly to remove any excess adhesive that may have been left behind. Apply the new part by removing the backing material and position as shown. Press firmly to ensure adhesion.

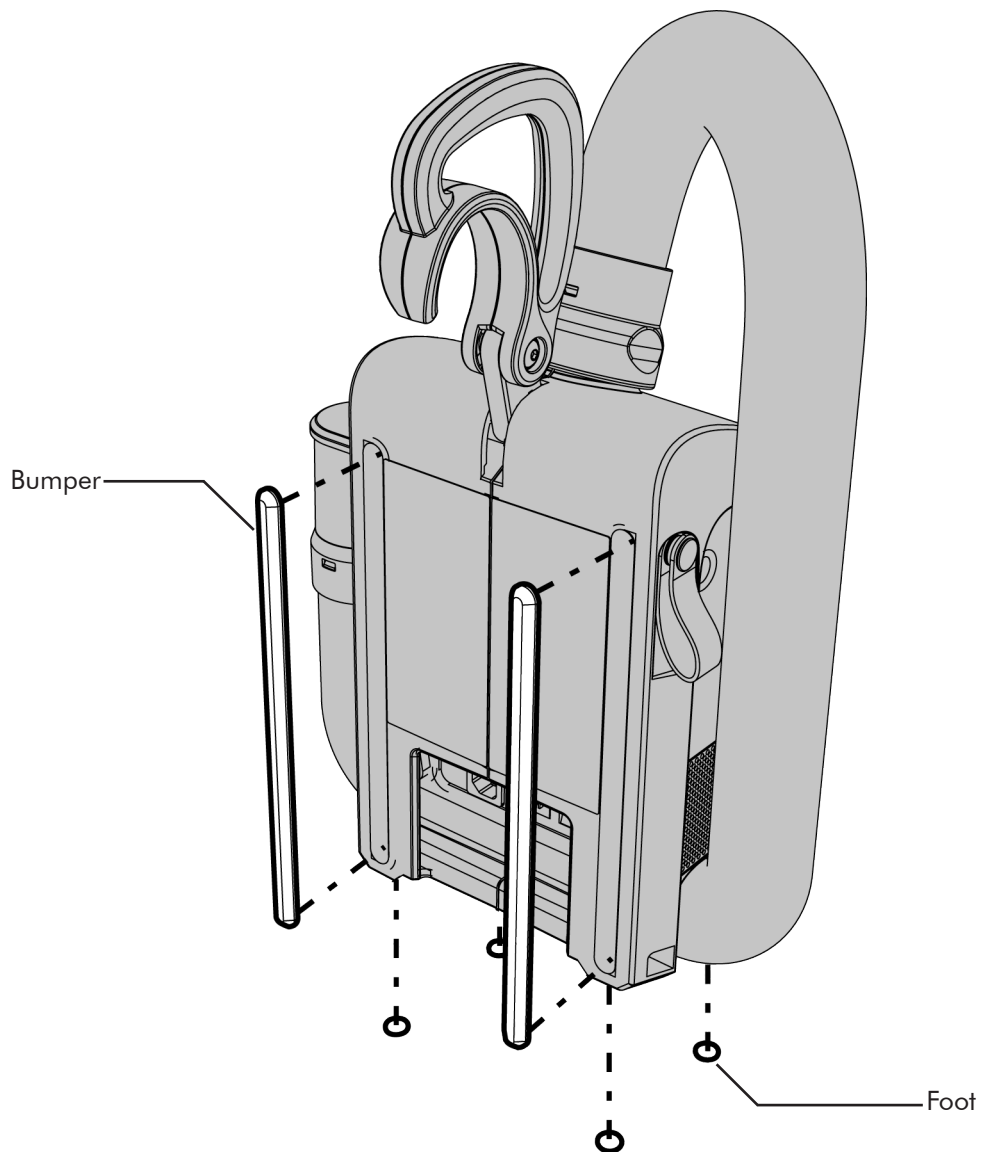

Hose Removal

1. Remove the damaged hose by lifting the dust cover slightly and unsnapping the hose from the side of the unit as shown.
2. Push the release button at the top of the unit to remove the hose.
3. Attach the new hose by reversing the process.

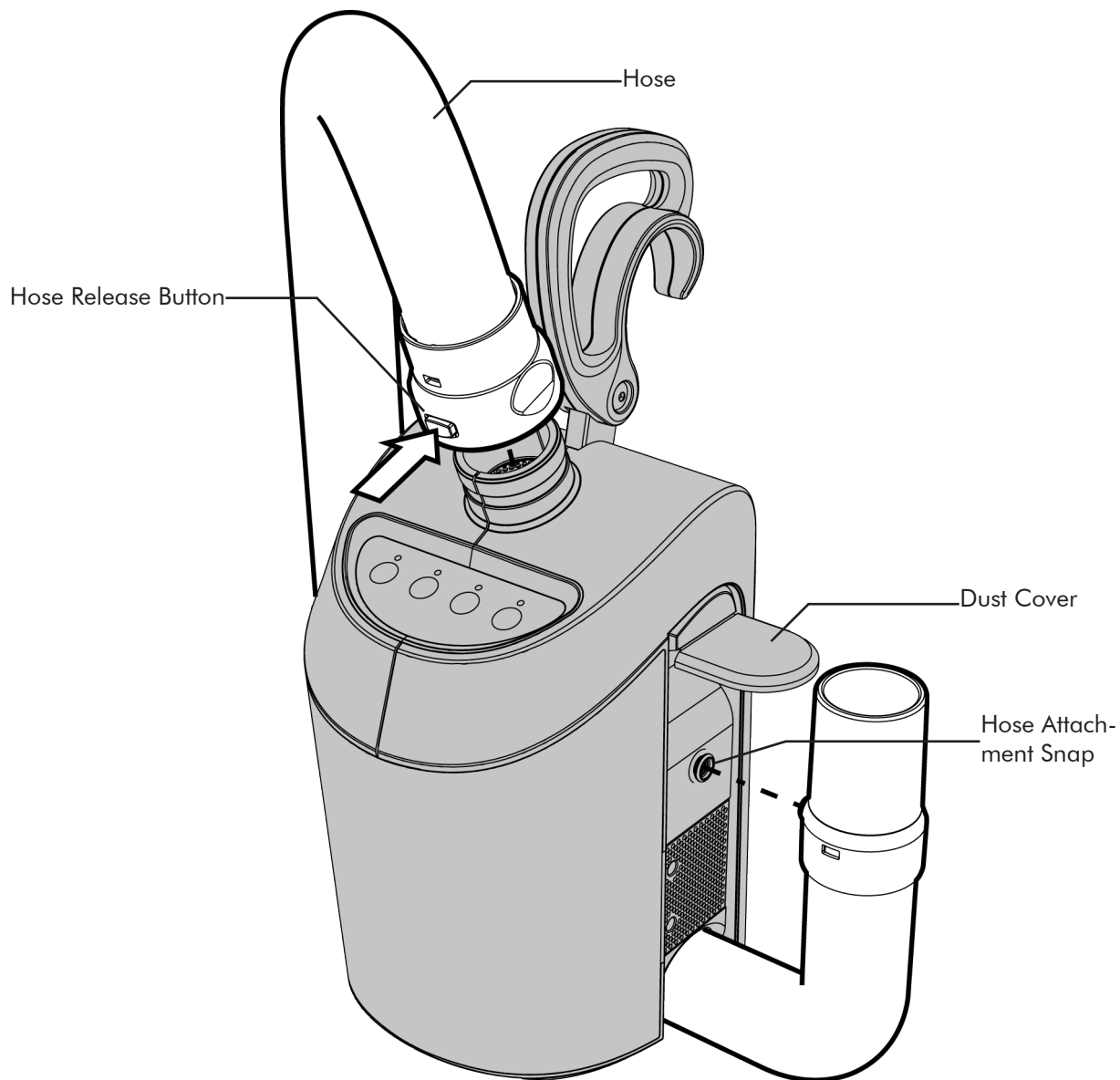

Air Filter and Air Filter Cover Replacement


To Remove Air Filter Covers

1. Disconnect hose from unit. (See page 23)
2. Remove the two phillips head screws on each side to detach the air filter covers.
3. Re-attach the new air filter covers and screws.

To Remove Air Filters

1. Follow the steps shown above to remove air filter covers.
2. Remove air filters and replace.
3. Re-attach the air filter covers and screws.

* For air filter cleaning instructions please refer to Page 16.

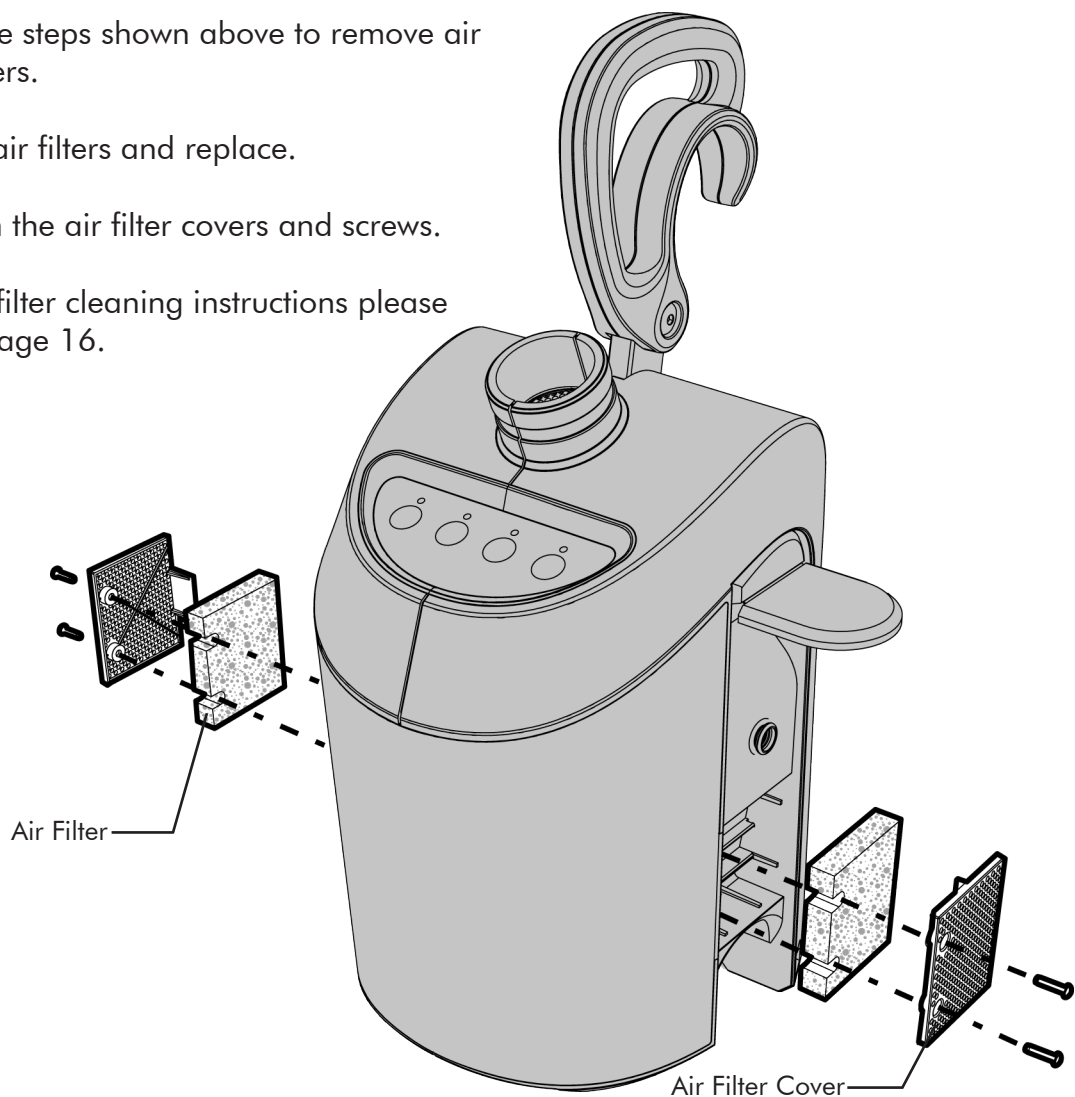


Dust Cover/Hose Attachment Snap Replacement

To Remove Dust Cover

1. Disconnect hose from unit. (See page 23)
2. Lift the cover "flap" to remove the 3 phillips head screws that attach the dust cover.
3. To replace dust cover, fold in the three flexible tabs. Then insert the screws one at a time starting with the top middle followed by the sides.


To Remove Hose Attachment Snap

1. Disconnect hose from unit. (See page 23)
2. Remove the phillips head screw and snap.

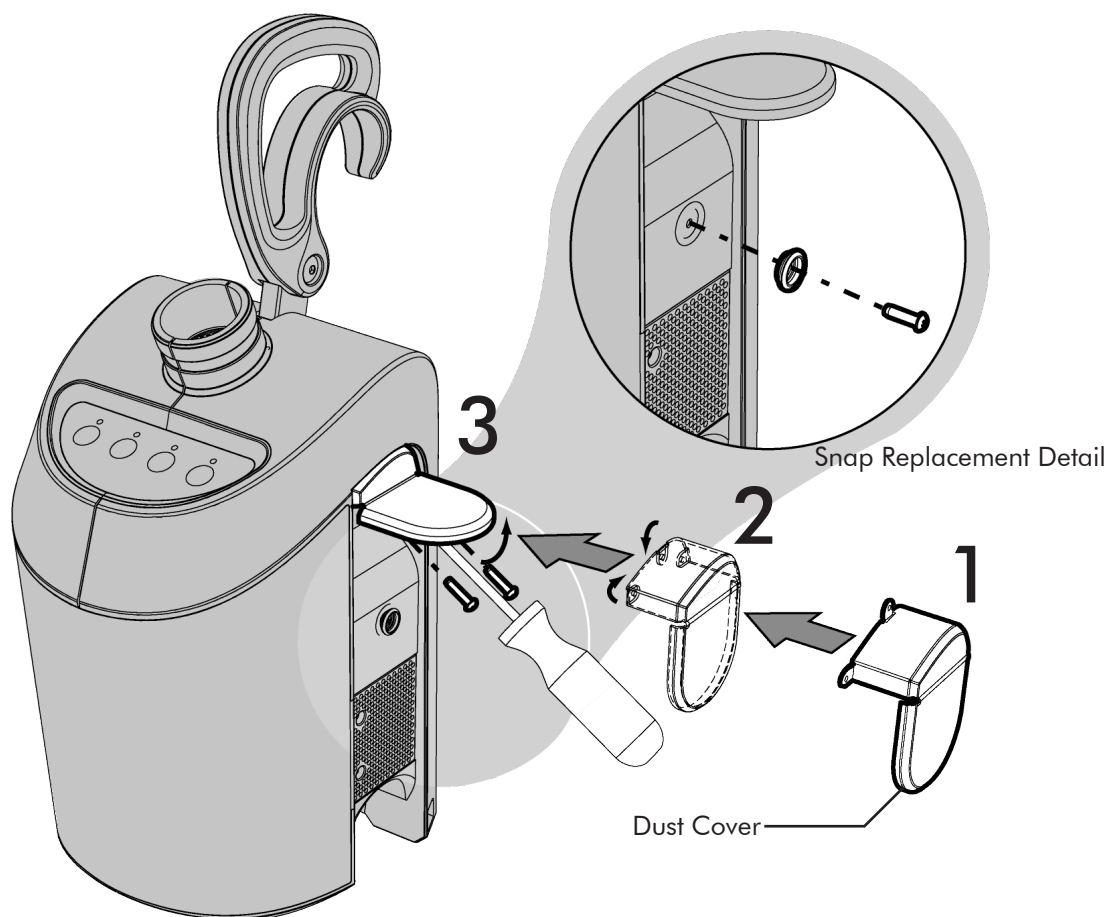


Metal Cover Replacement

1. Disconnect hose from unit. (See page 23)
2. Remove the two phillips head screws on each side to detach the air filter covers. Remove the air filters.
3. Gently pull the damaged cover apart to allow clearance to remove it.
4. Carefully re-insert the new metal cover. Reassemble the two filters and filter covers.
5. Secure everything in place by reinserting four screws.

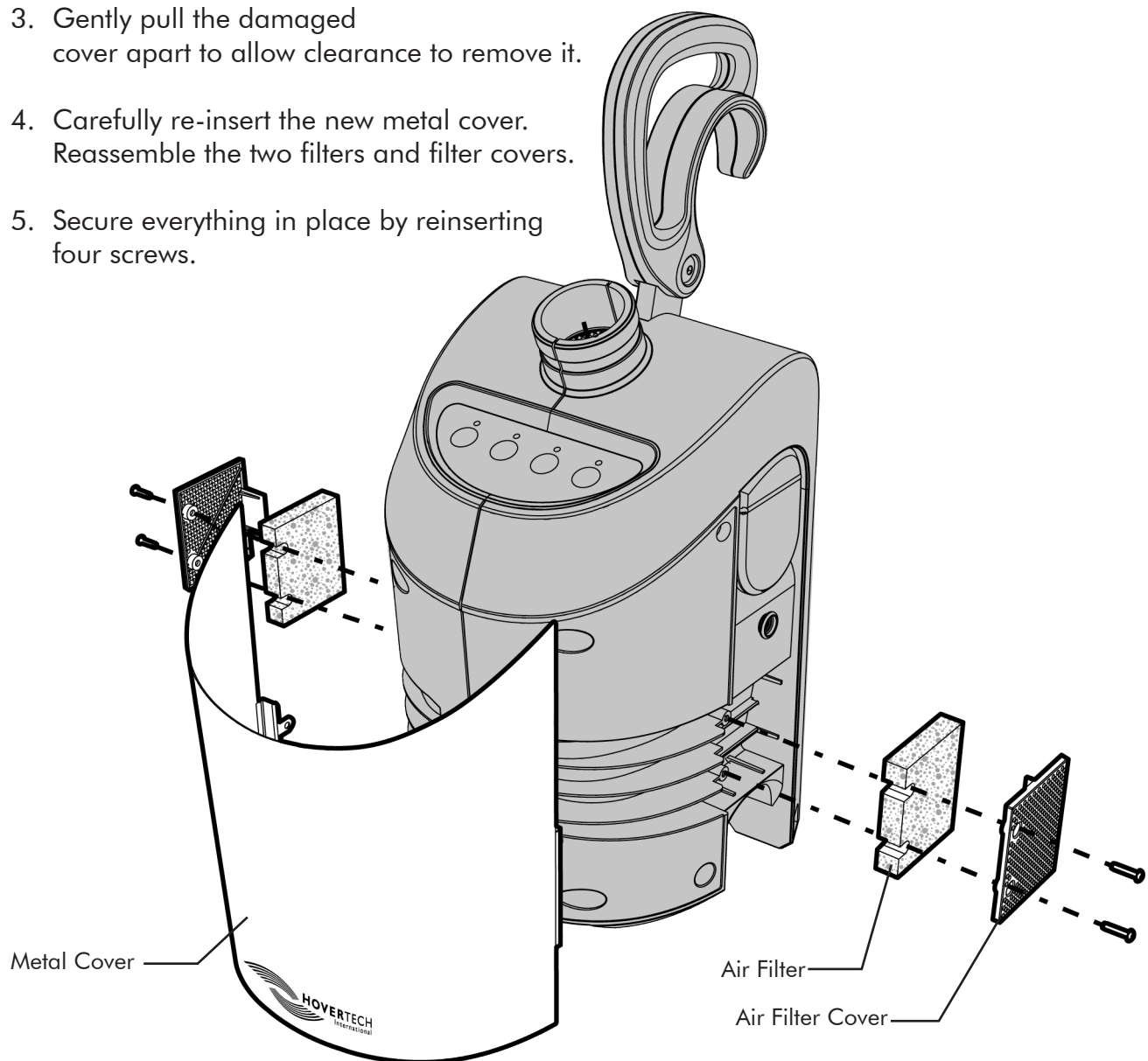

Cord Strap Replacement

1. Unsnap the strap and remove power cord.
2. Detach the damaged cord strap by removing the screw as shown.
3. Reattach strap by positioning in place and securing it with the screw provided.

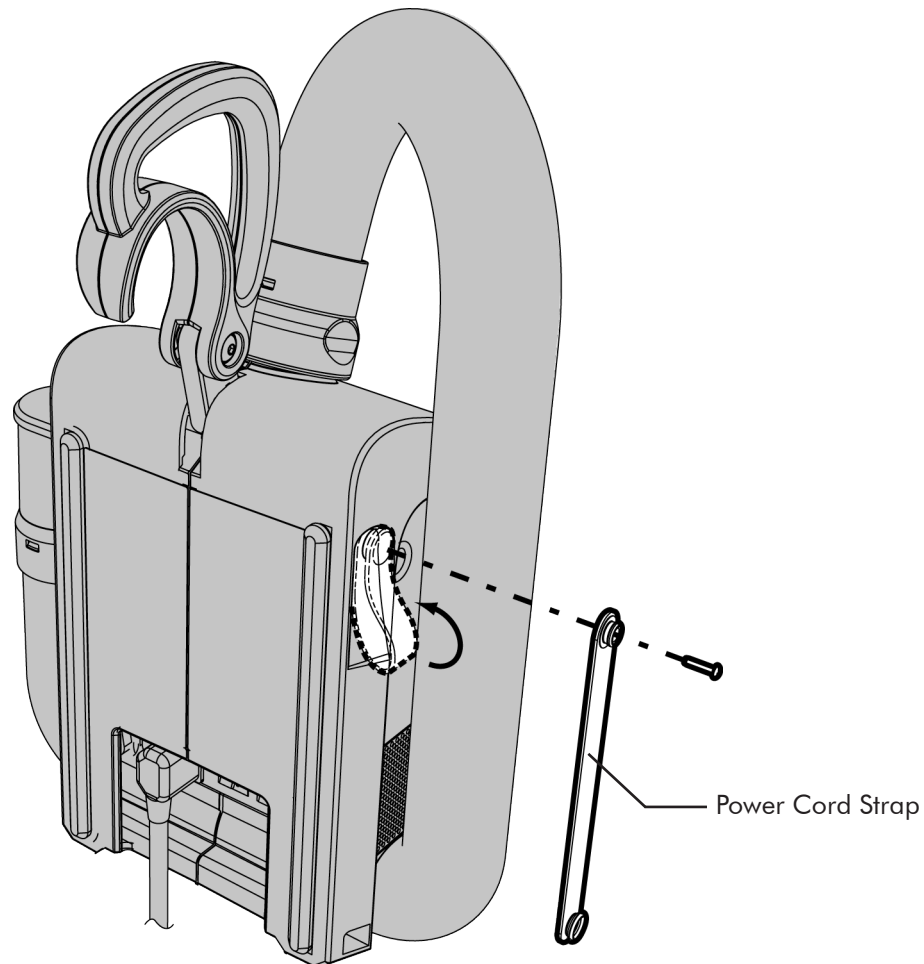

Troubleshooting

<div> <div>No Indicators (lights)</div> <div>No Air Flow</div> <div>No Button Function</div> <div>Mattress Pressure Low</div> </div>				
X	X	X		Verify AC from wall Check power cord connections at air supply and wall Check circuit breakers on rear of air supply unit
X				Return unit for repair
	X			Return unit for repair
		X		Return unit for repair
			X	Check hose connections at air supply and mattress Check hose for rips/tears Check mattress for rips/tears Verify that air filters are clean

HOVERMATT® USER MANUAL

Component Parts List

HoverTech Part Number	Description
HTA-AF	Air Filters (sold in pairs)
HTA-AFC	Air Filter Cover with screws (sold in pairs)
HTA-B/F	Bumpers and Feet (sold as a kit)
HTA-CPL	Control Panel Label
HTA-DC	Dust Cover w/ screws (kit)
HTA-H/AH	Handle/Attachment Hook
HTA-HAS	Hose Attachment Snap
HTA-Hose	Hose Assembly
HTA-IL	Information Labels
HTA-MC	Metal Cover w/ screws
HTA-PCord	Power Cord
HTA-PCS	Power Cord Strap with screw
HTA-PCC	Power Cord Clamp

Warranty Statement

The HoverMatt® air transfer system and the HoverTech International Air Supplies are warranted to be free from defects in materials and workmanship for (1) one year. Warranty begins on date of in-service by a HoverTech International representative or shipment date.

In the unlikely event that a problem arises as a result of a defect in materials or workmanship, we will promptly repair your item or replace it if we feel that it cannot be repaired – at our expense and discretion using current models or parts performing the equivalent function – upon receipt of the original item to our repair department. You must pre-notify HoverTech International by phone (800-471-2776). Item is to be insured by you against loss during transportation and must be shipped with transportation and/or broker charges prepaid. Should any HoverTech International product be returned that is not covered under warranty, there will be a minimum \$100 service charge plus shipping costs. Lead time for repairs is approximately 2 weeks. Please refer to the Return and Repairs section of this Manual for return instructions.

This warranty is not an unconditional guarantee for the life of the product. Our warranty does not cover product damage that may result from use contrary to Manufacturer's instructions or specifications, misuse, abuse, tampering, or damage due to mishandling. Warranty specifically does not cover product damage that may result from using an air supply that produces more than 3.5 psi to inflate the HoverMatt® air transfer system. Equipment that has been neglected, improperly maintained, repaired or altered by someone other than an authorized representative of Manufacturer, or operated in anyway contrary to the operating instructions, shall void this warranty.

This warranty does not cover normal "wear and tear". Component parts, particularly any optional equipment, valve caps, their attachments and cords, will show wear with use over time and eventually may need to be refurbished or replaced. This normal type of wear is not covered by our warranty, but we will provide prompt, high quality repair service and parts at a nominal cost.

Warranty Statement

HoverTech International's liability under this warranty and on any claim of any kind for any loss or damage arising out of, connected with, or resulting from the design, manufacture, sale, delivery, installation, repair or operation of its products, whether in contract or tort, including negligence, shall not exceed the purchase price paid for the product and upon expiration of the applicable warranty period, all such liability terminates. The remedies which this warranty provides are exclusive and HoverTech International shall not be liable for any incidental or consequential damages.

There are no warranties, expressed or implied, which extend beyond this warranty statement. The provisions of these warranty clauses are in lieu of all other warranties, expressed or implied, and of all other obligations or liabilities on HoverTech International's part and neither assumes nor authorizes any other person to assume for HoverTech International any other liability in connection with Manufacturer sale or lease of said products. HoverTech International makes no warranty of merchantability or fitness for a particular purpose. There is no warranty that the goods will be fit for a particular purpose. By accepting the goods, the buyer acknowledges that buyer has determined the goods are suitable for the buyer's purposes.

MANUFACTURER'S SPECIFICATIONS ARE SUBJECT TO CHANGE.

Returns and Repairs

All products being returned to HoverTech International must have a Return Goods Authorization Number issued from the company. Please call 800-471-2776 for an RGA #. Any products returned without the necessary RGA # may cause a delay in the repair time. If the product is not covered under warranty, a minimum charge of \$100 will be assessed for each repair. Should a repair charge be assessed, HoverTech International will notify the facility and a purchase order for the repair will need to be issued before the repair can be completed. Lead-time for repairs is approximately 2 weeks. A fee of \$50 will be assessed if HoverMatt® air transfer systems are not properly cleaned for repair.

All products should be sent to:

HoverTech International
513 South Clewell St.
Bethlehem, PA 18015

Attn: Repair Dept./RGA # _____
Phone: 800-471-2776
Fax: 610-694-9601



T: 1300 106 003

HoverTech International

513 South Clewell St.
Bethlehem, PA 18015
Phone: 800-471-2776
www.HoverMatt.com